

CEF-MHC Class I Control Peptide Pool "Classic"

Catalog #: CTL-CEF-001

Lot: _____

Product: CEF-MHC Class I Control Peptide Pool "Classic"

Formulation: Each vial contains a total of 460µg lyophilized peptide powder consisting of 23 individual peptides, 20µg of each peptide.

Description: The CEF-MHC Class I Control Peptide Pool "Classic" contains 23 peptides, each corresponding to a defined HLA class I-restricted T cell epitope from cytomegalovirus, Epstein-Barr virus and influenza virus.

Activity: The CEF-MHC Class I Control Peptide Pool "Classic" stimulates corresponding peptide-specific CD8+ memory T cells to release IFN-γ. These peptides are restricted by HLA-A1, -A2, -A3, -A11, -A24, -A68, -B7, -B8, -B27, -B35, and -B44 molecules, and have been shown to elicit recall responses in the majority individuals expressing these rather common MHC class I alleles – most humans have been previously exposed to these pathogens. Because of promiscuous binding motifs, these peptides can occasionally bind other MHC class I, but not class II molecules, further extending the range of individuals in whom they elicit a CD8 cell recall response. Therefore, the majority of randomly selected human donors will respond to the CEF-MHC Class I Control Peptide Pool "Classic".


Recommended Use: The CEF-MHC Class I Control Peptide Pool "Classic" is recommended as the positive control for detecting antigen-specific CD8+ cells in human PBMC, for example, when performing cytokine assays for immune monitoring purposes. Such assays include IFN-γ measurements by ELISPOT and intracytoplasmic cytokine staining (ICS), for accurate frequency measurements of the cytokine producing CD8 cells, or cytokine bead arrays (CBA), and ELISA, for a semi-quantitative read out.

Instruction for Use: **Stock solution:** Flick tube to ensure all powder is at the bottom of the tube. Add 10µl tissue culture grade DMSO followed by 40µl of sterile double distilled water. Vortex briefly and watch sterile handling! Add 450µl of tissue culture grade PBS and vortex briefly. Flick tube to ensure that all liquid is at the bottom of the tube. The stock solution is ready for use. This stock solution (500µl, 20x) contains each peptide at 40µg/ml concentration. The stock can be stored at 4°C for one week. CTL recommends to aliquot and store the stock solution at -20°C to -80°C for long term storage.

Use peptide pool at 2µg/ml of final peptide concentration that is at 1:20 dilution of the stock. For example, prepare a 1:10 dilution of the stock in tissue culture medium as the working solution, and add it 1+1 (v/v) to the medium containing the PBMC. For ELISPOT assays we recommend plating 100µl of this 2x peptide working solution per well, directly into the ELISPOT plate, followed by



Product Data Sheet

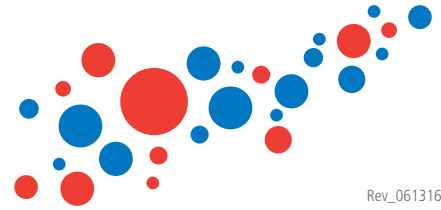

Rev_061316

adding the 100µl of the PBMC for a 24h assay duration. For exact frequency measurements, we recommend testing in triplicate. Because serum is the greatest variable in assay performance, we recommend the use of serum-free media at all steps of the assay. CTL offers such serum-free media (CTLT-005) that has been customized for low background/high signal performance with PBMC.

Storage: In lyophilized powder form at 4°C. The stock solution can be stored at 4°C for one week; however -20°C to -80°C is suitable for long-term storage.

Peptides: VSDGGPNLY (Influenza A, HLA-A1), CTELKLSY (Influenza A, HLA-A1), GILGFVFTL (Influenza A, HLA-A2), GLCTLVAML (EBV, HLA-A2), NLVPMVATV (HCMV, HLA-A2), KTGGPYKR (Influenza A, HLA-A68), ILRGSVAHK (Influenza A, HLA-A3), RVRAYTYSK (EBV, HLA-A3), RLRAEAQVK (EBV, HLA-A3), IVTDFSVIK (EBV, HLA-A11), ATIGTAMYK (EBV, HLA-A11), DYC NVLNKEF (EBV, HLA-A24), RPPIFIRRL (EBV, HLA-B7), ELRSRYWAI (Influenza A, HLA-B8), RAKFKQLL (EBV, HLA-B8), FLRGRAYGL (EBV, HLA-B8), QAKWRLQTL (EBV, HLA-B8), SRYWAIRTR (Influenza A, HLA-B27), RRIYDLIEL (EBV, HLA-B27), YPLHEQHGM (EBV, HLAB35), EENLLDFVRF (EBV, HLA-B44), EFFWDANDIY (HCMV, HLA-B44), TPRVTGGGAM (HCMV, HLA-B7).

References: Currier JR, Kuta EG, Turk E, Earhart LB, Loomis-Price L, Janetzki S, Ferrari G, Birx DL, Cox JH. A panel of MHC class I restricted viral peptides for use as a quality control for vaccine trial ELISPOT assays. *J Immunol Methods* 260:157-172, 2002.

Mwau M, Michael A and Hanke T. Design and validation of an enzyme-linked immunospot assay for use in clinical trials of candidate HIV vaccines. *AIDS Res Hum Retroviruses* 18: 611-8, 2002.

For laboratory research use only. Not for use in diagnostic or therapeutic procedures.